






Engage in the future

### Discover the benefits of being connected

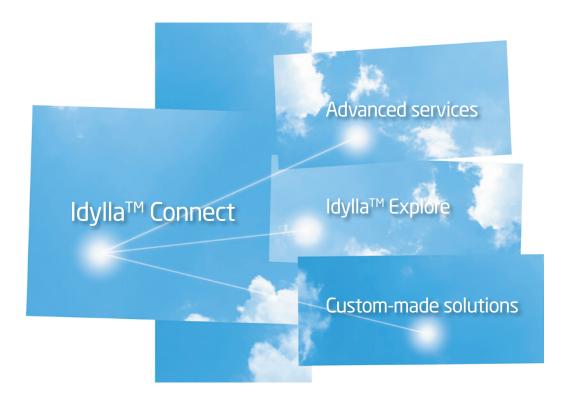

#### Advanced services to ensure continuity in your laboratory workflow



# Automatic software updates

New releases of assay and console software are sent to your  $dylla^{TM}$  console and can be installed with a single touch on the screen.



# Immediate and remote service and support

ldylla<sup>™</sup> system parameters and error logs can be analyzed at anytime and anywhere to ensure swift actions and solutions.

# What about privacy and data protection?

 $(\begin{subarray}{c} \end{subarray})$ 

Is Idylla™ Connect safe?

The Idylla<sup>™</sup> Connect network is built on state of **the art cyber security technology** which guarantees that only Idylla<sup>™</sup> Systems can get access to this network. Connecting to Idylla<sup>™</sup> Connect requires a **certificate based authentication procedure** which ensures that only Biocartis known and trusted Systems can connect. The Idylla<sup>™</sup> Connect and Idylla<sup>™</sup> Explore solution complies with the US FDA adopted cyber security standards of the CLSI organization, AUTO-9A (Remote Access to Clinical Laboratory Diagnostic Devices via the Internet) and AUTO11-A2 (Remote Access to Clinical Laboratory Diagnostic Devices via the Internet).



Is my data secure?

Data is transported over a single TLS (Transport Layer Security) socket between your Idylla<sup>TM</sup> System and the Biocartis server to provide **secure end-to-end confidentiality of encrypted data.** The Idylla<sup>TM</sup> Connect data transfer technology prevents tampering with the data.



What about Patient Privacy?

Idylla<sup>™</sup> does not require patient identifying information. Idylla<sup>™</sup> test results only reference to an **anonymous sample ID**. Biocartis recommends users to avoid including patient names or other medical information in the sample ID or annotations field.



Who can access my data?

Access to data on Idylla<sup>TM</sup> Explore is protected with **login credentials** that are associated to specific data group definitions. This architecture ensures that test result data is **only accessible by authorized persons** within the group. The Biocartis server, holding the transferred test result data, is physically located at the European Biocartis headquarters and maintained by Biocartis staff members with audit functionality.




#### More insight into your data with Idylla<sup>™</sup> Explore

Get connected and enjoy the advantages of Idylla<sup>™</sup> Explore, a web-based application that allows you to analyze your data by providing

- Visualization of PCR curves from Idylla™Test Results
- Cq values per target
- Direct Access to Console result reports

Idylla<sup>™</sup> Explore can be accessed anywhere and anytime from your PC or tablet through the following link: https://idyllaexplore.biocartis.com

Subscribe today and join the Idylla<sup>™</sup> Explore community by sending an email to explore@biocartis.com

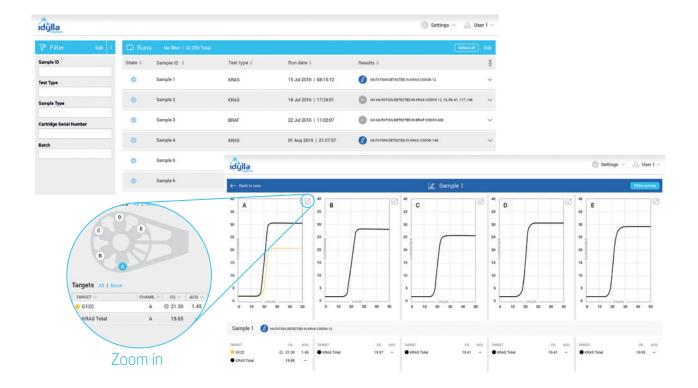

### Connect today

Get started with Idylla<sup>™</sup> Connect and Idylla<sup>™</sup> Explore in just 2 easy steps




Establish an ethernet connection between the ldylla<sup>™</sup> console and your laboratory network or a 3/4G modem



On the Idylla<sup>TM</sup> console:

- Enable Idylla<sup>™</sup> Connect
- Enable "Allow uploading of test result data Idylla™ Explore" for Idylla™ Explore

## Your data is safe with us

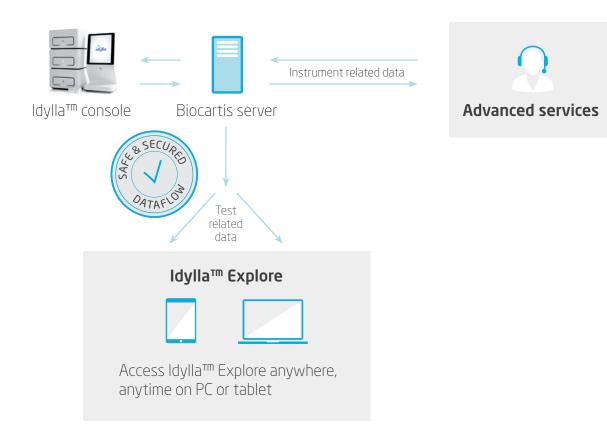

#### Custom-made solutions

Biocartis offers personalised solutions that fit your individual needs

- Create a network between different Idylla™ User sites and share data and knowledge
- Direct access to your data for building your own solution
- Statistical analysis on your obtained data
- Monitoring: Follow-up of your data over time
- Disease surveillance or diagnostic grid: linking of real-time molecular diagnostic test data to geo-location and sample data

# Would you like to know more?

Your local sales representative will be more than happy to inform you.

#### Contact

Biocartis NV Generaal De Wittelaan 11 B3 2800 Mechelen – Belgium T +32 15 632 800

#### customersupport@biocartis.com www.biocartis.com

Biocartis trademark and logo are trademarks belonging to Biocartis and are used and registered in Europe. Idylla is a registered trademark in the United States and other countries. Idylla trademark and logo are used trademarks belonging to Biocartis.

CONFIDENTIALITY STATEMENT: Biocartis will not require access to confidential patient-identifying information. Biocartis assumes that the sample ID field does not contain any patient-identifying information. Should that be the case, Biocartis will treat it as confidential information and will not use it for any purpose.

ldylla<sup>™</sup> Explore has been developed solely as a visualization tool and does not allow any modification of reported data. Idylla<sup>™</sup> Explore may not be available in all markets due to regulatory restrictions. Specific terms and conditions apply for the use of Idylla<sup>™</sup> Explore.





